




Temperature Module GEN2

Quickstart Guide

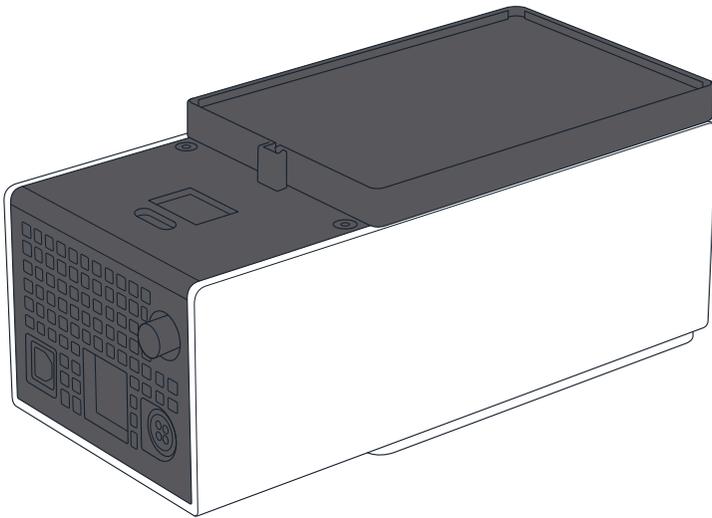

Opentrons Labworks Inc.

October 2023


Product Description

The Opentrons Temperature Module GEN2 is a hot and cold plate module. It can reach and maintain temperatures ranging from 4 °C to 95 °C. The module is compatible with the Opentrons Flex and Opentrons OT-2 liquid handling robots.


Table of Contents

Product Elements	3
– Box Contents	
– Thermal Blocks	
– Flex Caddies	
Before You Begin	6
– Deck Placement and Cable Alignment	
– Anchor Adjustments	
Flex Attachment Steps	9
OT-2 Attachment Steps	10
Additional Product Information	10
– Maintenance	
– Warranty	
– Support	
– App Download	
– Completed Certifications	
– Recommended Operating Conditions	
– Manufacturer Description	

Product Elements

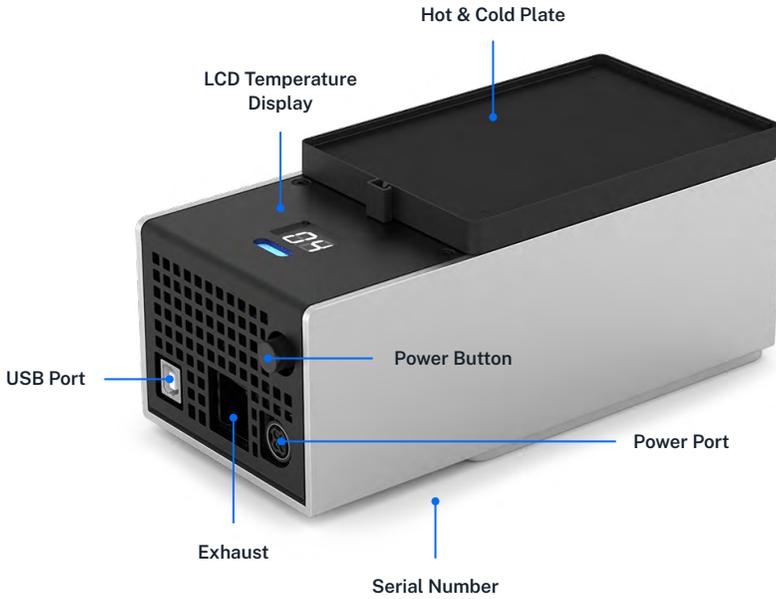


BOX CONTENTS

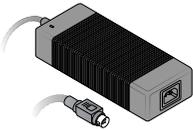

(1) Power Supply



(1) Power Cable

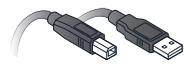


(1) USB Cable



(1) 24-Well Thermal Block



(1) 96-Well PCR Thermal Block

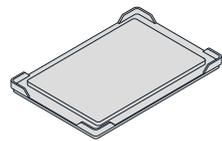


(1) Flat Bottom Plate for OT-2


THERMAL BLOCKS

The Temperature Module comes with 24-well, 96-well, and flat bottom thermal blocks. The blocks hold 1.5 mL and 2.0 mL tubes, 96-well PCR plates, PCR strips, and flat bottom plates.

The Temperature Module Caddy comes with a deep well thermal block and a flat bottom thermal block designed for use with the Flex Gripper. These blocks hold deep well plates and flat bottom plates.

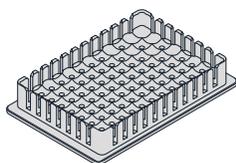

Deep Well Thermal Block

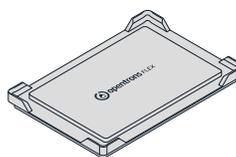

Flat Bottom Plate for Flex


Opentrons recommends only using the following thermal blocks, depending on the type of robot you have.

Thermal block	Flex	OT-2
24-well	✓	✓
96-well PCR	✓	✓
Deep well	✓	✗
Flat bottom for Flex	✓	✗
Flat bottom for OT-2	✗	✓

You can tell which flat bottom thermal block you have because the one for Flex has the words “Opentrons Flex” on its top surface. The one for OT-2 does not.

FLEX CADDIES

When used with a Flex robot, the Temperature Module GEN2 fits into a caddy that occupies space below the deck. The caddy places your labware closer to the deck surface and allows for below-deck cable routing. See the Modules chapter in the [Flex Instruction Manual](#) for more information.

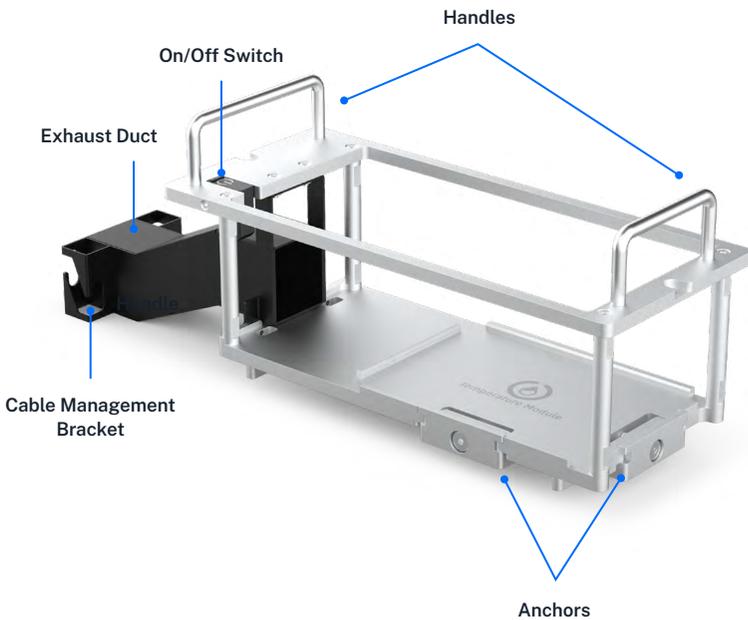

The OT-2 does not use caddies. Modules clip directly to the deck.

Module caddies are available for purchase at shop.opentrons.com.


Before You Begin

Review this section for important information about Temperature Module GEN2 deck placement, alignment, and anchor adjustments before installing the module.

DECK PLACEMENT AND CABLE ALIGNMENT

The supported deck slot positions for the Temperature Module GEN2 depend on the robot you're using.

Robot model	Deck placement
Flex	In any deck slot in column 1 or 3. The module can go in slot A3, but you need to move the trash bin first.
OT-2	In deck slots 1, 3, 4, 6, 7, 9, or 10.

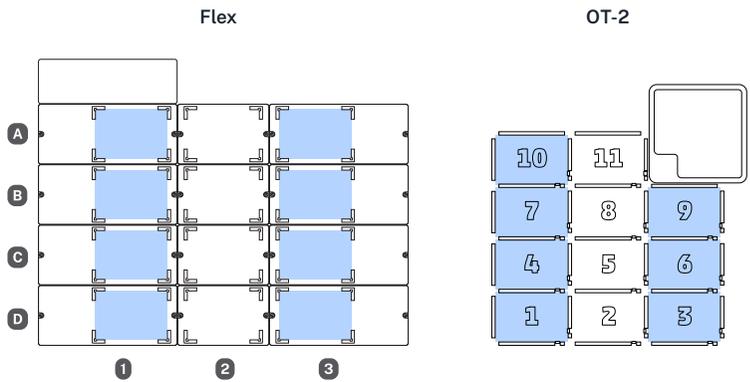

To properly align the module relative to the robot, make sure its exhaust, power, and USB ports face outward, away from the center of the deck. This keeps the exhaust port clear and helps make cable routing and access easier.

Robot model	Exhaust, Power, and USB alignment
Flex	<ul style="list-style-type: none">Facing left in column 1.Facing right in column 3.
OT-2	<ul style="list-style-type: none">Facing left in slot 1, 4, 7, or 10.Facing right in slot 3, 6, or 9.




Warning: Do not install the Temperature Module GEN2 with the ports facing in, towards the middle of the deck. This alignment vents air onto the deck and makes cable routing and access difficult.

ANCHOR ADJUSTMENTS

Anchors are screw-adjustable panels on the Temperature Module GEN2 caddy. They provide the clamping force that secures the module to its caddy. Use the 2.5 mm screwdriver to adjust the anchors.

- To loosen/extend the anchors, turn the screws counter-clockwise.
- To tighten/retract the anchors, turn the screws clockwise.

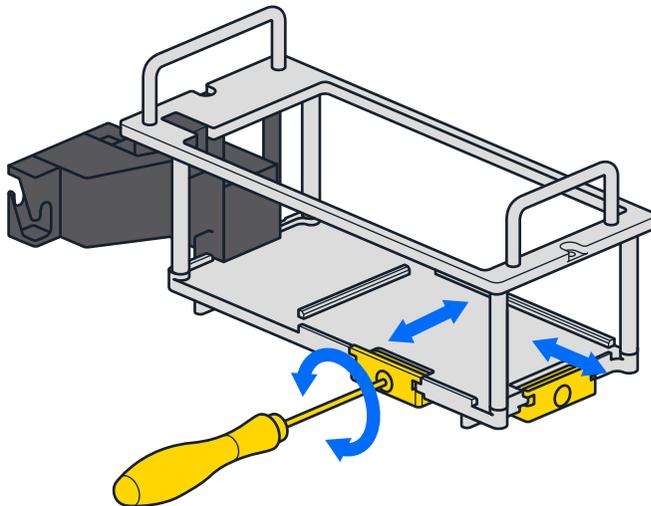

Before installation:

- Check the anchors to make sure they're level or extend slightly past the base of the caddy.
- If the anchors interfere with installing the module, adjust them until there's enough clearance to seat the module and then tighten them to hold it in place.

Flex Attachment Steps

To attach the Temperature Module GEN2 to your Flex:

1. Choose the supported slot you want to use for the module. Use the 2.5 mm screwdriver to remove the deck slot plate.
2. Insert the module into its caddy by aligning the power button on the module with the on/off switch on the caddy.



Tip: If you're having trouble inserting the module into its caddy, the module's power button is probably facing away from the caddy's on/off switch. Turn the module so the power button faces the on/off switch and try again.

3. Holding the module in the caddy, use the 2.5 mm screwdriver to turn the anchor screws clockwise to tighten the anchors. The module is secure when it doesn't move while gently pulling on it and rocking it from side to side.
4. Connect the power and USB cables to the module and route the cables through the cable management bracket on the end of the caddy's exhaust duct.
5. Insert the caddy into the deck slot, exhaust duct first, and route the power and USB cables through the Flex. *Do not* connect the power cable to a wall outlet yet.
6. Connect the free end of the USB cable to a USB port on the Flex.
7. Connect the power cable to a wall outlet.
8. Gently press the on/off switch to turn the module on. If you successfully connected the module, it will appear in the Pipettes and Modules section on your robot's device detail page in the Opentrons App.

OT-2 Attachment Steps

To attach the Temperature Module GEN2 to your OT-2:

1. Choose the supported slot you want to use for the module and press it gently into place.
2. Connect the USB cable to the module and to a USB port on the OT-2. Use the cable holders on the side of the robot to keep the USB cable out of the way.
3. Connect the power cable to the module and then to a wall outlet.
4. Gently press the on/off switch to turn the module on. If you successfully connected the module, it will appear in the Pipettes and Modules section on your robot's device detail page in the Opentrons App.

Additional Product Information

MAINTENANCE

Users should not attempt to service or repair the module themselves. If you have concerns about your module's performance or require maintenance, please contact Opentrons Support.

WARRANTY

All hardware purchased from Opentrons is covered under a 1-year standard warranty. Opentrons warrants to the end-user of the products that they will be free of manufacturing defects due to part quality issues or poor workmanship and also warrants that the products will materially conform to Opentrons' published specifications.

SUPPORT

Opentrons Support can help you with questions about our products and services. If you discover a defect, or believe your product is not functioning to published specifications, contact us at support@opentrons.com.

Please have the Temperature Module's serial number available when contacting support. You can find the serial number on the bottom of the module or in the Opentrons App. On the Temperature Module card in the Pipettes and Modules section of your robot's device detail page, click the three-dot menu (:) and then **About**.

APP DOWNLOAD

Control your liquid handling robot and module using the Opentrons App. Download the app for Windows, macOS, or Ubuntu at <https://opentrons.com/ot-app/>.

COMPLETED CERTIFICATIONS

IEC, FCC, RoHS

RECOMMENDED OPERATING CONDITIONS

Environmental temperature: 20–25 °C
Environmental humidity: 80% maximum

MANUFACTURER DESCRIPTION

Opentrons Labworks Inc
45-18 Ct Square W
Long Island City, NY 11101



Post-sales service & contacting Opentrons

If you have any questions about the use of the system, abnormal phenomena, or special needs, please contact: support@opentrons.com. Also visit www.opentrons.com.